velleman®

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INFRARED THERMOMETER
INFRAROOD THERMOMETER
THERMOMÈTRE INFRAROUGE
TERMÓMETRO INFRARROJOS
IR-THERMOMETER

USER MANUAL
GEBRUIKERSHANDLEIDING
NOTICE D'EMPLOI
MANUAL DEL USUARIO
BEDIENUNGSANLEITUNG

CE

SHE11 – INFRARED THERMOMETER

1. Introduction

To all residents of the European Union

Important environmental information about this product

This symbol on the device or the package indicates that disposal of the device after its lifecycle could harm the environment. Do not dispose of the unit (or batteries) as unsorted municipal waste; it should be taken to a specialized company for recycling. This device should be returned to your distributor or to a local recycling service.

Respect the local environmental rules.

If in doubt, contact your local waste disposal authorities.

Thank you for choosing Velleman! Please read the manual thoroughly before bringing this device into service. If the device was damaged in transit, do not install or use it and contact your dealer.


The **SHE11** is an easy-to-use ear thermometer with high-speed temperature reading, ambient temperature reading, alarm, clock and memory function.

2. Safety Instructions

- Damage caused by disregard of certain guidelines in this manual is not covered by the warranty and the dealer will not accept responsibility for any ensuing defects or problems.
- Clean the probe with a clean tissue and some rubbing alcohol. Do not immerse the thermometer in liquids.
- Do not switch the device on immediately after it has been exposed to changes in temperature. Protect the device against extreme temperatures, shocks and abuse.
- For accurate measurements, do not measure the body temperature within 30 minutes after meals, exercise or a bath.
- Note that damage caused by user modifications to the device is not covered by the warranty.

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3. Description

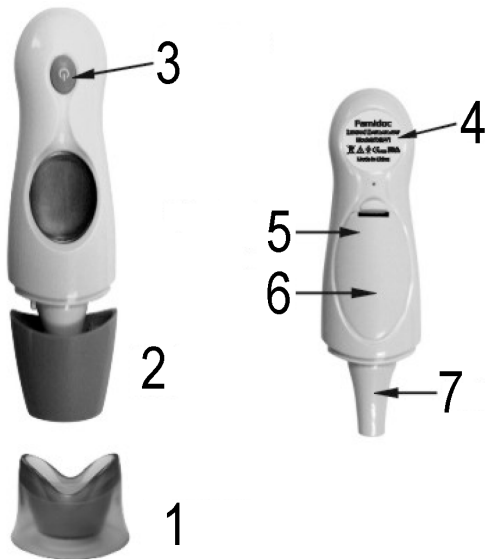- 1. base
- 2. cap
- 3. power/scan button
- 4. safety label

- 5. time setting button
- 6. battery cover
- 7. probe

4. Description of the Symbols

Display

Measuring in progress

Ear scan mode

Forehead scan mode

Fever

Celsius/Fahrenheit scale

Clock

Low battery

Last memorized temperature display

5. The LCD

Display	Meaning	Action
	Alternate time and temperature display at power-on.	Set time at first use or after installing battery.
	Initialization process.	Wait for the last memorized temperature to display.
	Last memorized temperature is displayed.	Thermometer is ready for measurement.
	Battery is running low.	Replace the battery.
	Insufficient power.	Replace the battery immediately.
	Measured temperature is too high.	Measured temperature is higher than 42.9°C or 109.2°F
	Measured temperature is too low.	Measured temperature is lower than 32°C or 89.6°F.
	Environment temperature exceeds working temperature scope.	Move the thermometer to a room at ambient temperature.

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6. How to Use the Thermometer

- **Time and Ambient Temperature**

Set the time at first use or after installing a new battery. Open the battery compartment, press the time setting button. Set the hours with the scan button and confirm with the time setting button. Now, set the minutes with the scan button and confirm with the time setting button.

When the time is set and confirmed, the thermometer will alternatively display the time and ambient temperature.

- **Taking Someone's Temperature – In-Ear**

Press the scan button. The last measured temperature will be displayed. Wait for the thermometer to beep twice. Now, remove the cap. Insert the thermometer probe in the ear and press the scan button. Wait for the beep and remove the thermometer. Read the measured value from the display.

A measured temperature exceeding 38.0°C or 100.4°F will be indicated by one long and two short beeps, and the fever icon.

- **Taking Someone's Temperature – Forehead**

Press the scan button. The last measured temperature will be displayed. Wait for the thermometer to beep twice. Leave the cap onto the thermometer. Place the thermometer against the forehead and press the scan button. Wait for the beep and remove the thermometer. Read the measured value from the display.

A measured temperature exceeding 38.0°C or 100.4°F will be indicated by one long and two short beeps, and the fever icon.

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• Notes

Clean the thermometer probe after use and replace the cap on the thermometer.

The temperature taken on someone's forehead is for reference only since the environment, sweat, etc. can heavily influence the reading.

Switch on the thermometer. Press the scan button and wait 3 seconds to select the °C or °F measuring unit.

• Battery Replacement

For accurate measuring, replace the battery as soon as the battery symbol is displayed. Open the battery compartment, take out the old battery and replace with 1 x CR2032 lithium battery. Insert the battery with the positive pole facing upwards.

7. Technical Specifications

Measuring Range	0.0°C ~ 50.0°C or 32.0°F ~ 122.0°F
Display Resolution	0.1°C or 0.1°F
Accuracy	± 0.2°C or ± 0.4°F for 36.0°C ~ 39.0°C or 69.8°F ~ 102.2°F range ± 0.3°C or ± 0.5°F for all other ranges
Measuring Unit	°C or °F
Operating Temperature	10°C ~ 40°C or 50°F ~ 104°F, ≤ 95% RH
Storing Temperature	-20°C ~ 50°C or -4°F ~ 122°F, ≤ 95% RH
Power Supply	1 x 3V CR2032C (order code CR2032C , incl.)
Weight	100g
Certifications	EN 12470-5:2003 ASTM E 1965-98 IEC/EN 60601-1 IEC/EN 60601-1-2 ISO 14971:2000 CE0123 FDA 510(K)052849

Use this device with original accessories only. Velleman nv cannot be held responsible in the event of damage or injury resulted from (incorrect) use of this device. For more info concerning this product, please visit our website www.velleman.eu. The information in this manual is subject to change without prior notice.

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